

COMPANION MINI MASSAGE GUN

INSTRUCTION MANUAL

EMA-C025

READ ALL INSTRUCTIONS CAREFULLY BEFORE USE.


IN THE BOX

- EMA-C025 Companion Mini Massage Gun (x1)
- Charging Cable (x1)

SAFETY PRECAUTIONS

- This is a non-medical product suitable for relieving muscle strain after exercising.
- Please use it according to your own needs and under the guidance of your doctor.
- This device is NOT a toy, keep away from children and pets.
- Only use the device according to the instructions in this manual. Never use the device for anything other than its intended purpose.
- Do not open this device or attempt to repair the unit yourself if it is not working properly. Doing so risks fire, shock or injury.
- When not in use, the device should be stored in a dry, cool location, away from sources of water or heat.
- Do not expose device or charging cable to extreme temperatures (heat or cold), open flames, humid or wet conditions.
- Do not submerge in water. Do not use in a location where the device or charging cable could fall into liquid.
- Unplug the charging cable when not in use. Do not pull the charging cable to disconnect, or use the cable as a handle.
- Do not use if device or charging cable is damaged. If device or charging cable is damaged, it MUST be serviced or replaced by the manufacturer or its service agent before using again.
- Low battery may cause dysfunction.
- Do not over charge battery.
- This device contains a built-in battery and should not be disposed of with regular household waste. Please consult the local authorities on the correct way to dispose of this device.

PRODUCT DIAGRAMS

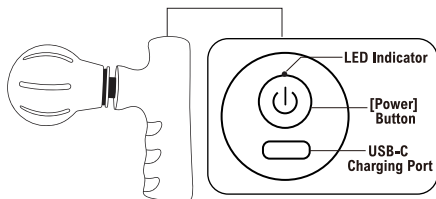

INSTRUCTIONS

Be sure to **fully** charge the speaker before using for the first time.

CHARGE THE MESSAGE GUN

- Plug the **small** connector of the included charging cable into the USB-C charging port on top of the massage gun.
- Plug the **large** connector into a USB power source, such as the USB port on a computer or into a compatible 5V USB adapter (not included) and plug into a wall socket.
- The LED indicator will **illuminate RED** when charging and turn OFF when fully charged.
- The massage gun **CAN NOT** be used during charging. It will power off automatically when being charged.

USING THE MESSAGE GUN

- Press and hold the **[Power] button** for 3 seconds to power on. The LED indicator will **glow BLUE** when massage gun is powered on.
- Press the **[Power] button** to change the toggle speed level between low, medium and high.
- Press and hold the power button for 3 seconds again to power off.
- The LED indicator will **flash BLUE** to warn of a low power of battery.

- The massage gun will power off automatically after 30 minutes on the same speed level.

NOTE:

- DO NOT use the massage gun directly on soft skin surface to avoid damage to your skin. Wear dry and soft clothing or place a dry towel on muscle area to be massaged.
- The massage time for the same area should not exceed 60 seconds.
- Press gently into the muscle. If not feeling pain or discomfort, you can apply a higher speed and add more pressure on the soft tissue parts of the body. DO NOT do this on the head or any bone parts of the body.
- Frequent high-speed and high-pressure massage on the same part will cause abrasions. Stop using it if you feel pain or discomfort.

TECH SPECIFICATIONS

- Battery Capacity: 1200mAh
- USB-C Input: DC 5V/1A
- Continuous Use Time:
 - Up to 1 hours - High Speed
 - Up to 2 hours - Medium Speed
 - Up to 3 hours - Low Speed
- Charging Time: around 2.5 hours
- Working Temperature: 14°F ~ 140°F (-10°C ~ 60°C)

FCC INFORMATION

This device complies with Part 15 of the FCC Rules.

Caution: Any changes or modifications not expressly approved could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio

frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.



Made in China